

BZA 100

Battery Impedance Analyser



- Impedance measurement of energy device ; battery , fuel cell or supercapacitors
- DC voltage measurement up to **100V**
- Quick diagnosis of test cells
- Battery lifetime estimation
- LAN interface with PC
- ZMAN impedance analysis software
- Cell temperature monitoring
- No requirement of additional electronic load or power source

Electrochemical impedance spectroscopy (EIS) is a widely used experimental technique to gain a deeper insight into the electrochemical processes of batteries. EIS cannot only provide detailed kinetic information, but can also be used to monitor changes in battery properties. EIS is a very sensitive technique, and offers a useful information about battery systems such as :

- Battery lifetime
- Internal temperature
- Internal defect

State of charge (SoC) represents the available battery capacity and is one of the most important states that need to be monitored to optimize the performance and extend the lifetime of batteries. Meanwhile, the State of Health (SoH) is an indicator associated with the long-term cycle life of batteries. EIS is increasingly being used for estimation of the SoH & SoC of batteries and it can be approximated with circles on the Nyquist plane, where the imaginary and the real part of the impedance are plotted on the Y-axis and X-axis, respectively. According to Fig. 1, R_s is quantified as the horizontal distance between the zero and the point where the EIS spectrum crosses the real axis (high frequencies) and R_s , R_{CT} and R_w are calculated as the horizontal distances of the each depressed semicircle, respectively. In Fig. 2, C_{DL} represent a double layer capacitance, R_{CT} is the charge transfer resistance, R_w stands for the Warburg impedance and R_s is the electrolyte resistance.

As shown in Fig. 3, the R_s value increases as performance deteriorates due to aging of the battery. Therefore, R_s is a criterion for determining the SoH of the battery. On the other hand, R_{CT} is related to the electrochemical reaction rate in the battery, and C_{DL} represents the double layer capacitance between the electrode and the electrolyte. As the depth of discharge (DoD) of the battery increases, SoC decreases. In the Nyquist plot, it can be seen that the diameter of the semicircle increases (R_{CT} increase. Fig. 4). Therefore, R_{CT} and C_{DL} are the criteria for indicating the SoC of the battery.

BZA 100 Battery Impedance Analyser

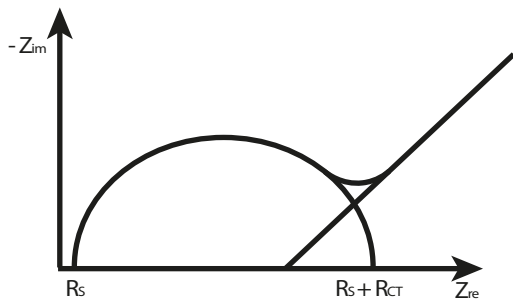


Fig. 1) Typical Nyquist plot

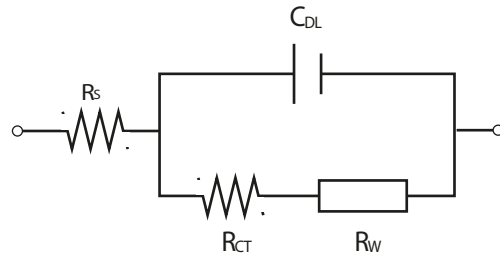


Fig. 2) Battery equivalent circuit model sample

C_{DL} : double layer capacitance

R_{CT} : charge transfer resistance

R_W : Warburg impedance

R_s : electrolyte resistance

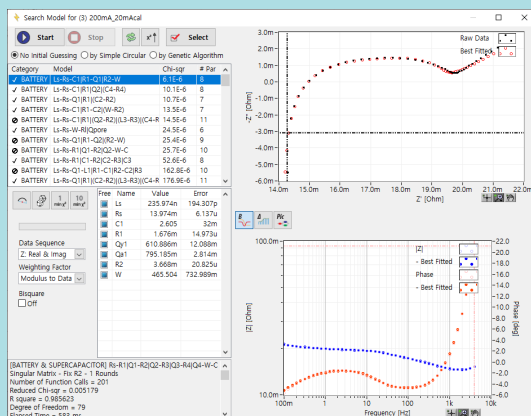


Fig. 3) Automatic Model Searching with ZMAN
BZA100 impedance data using 18650 battery

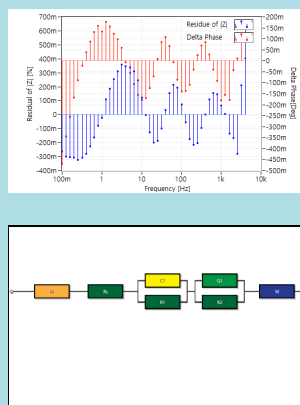


Fig. 4) 4kHz to 50mHz EIS
20points per decade

The **BZA100** Battery Impedance Analyzer, which covers a broad range of battery test functions ranging from DC voltage (up to 100V) and impedance test (500uΩ ~ 50Ω), is an ideal test tool for performance testing of single cell, battery pack ,individual stationary batteries under 100Volt.

The **BZA100** is designed to measure battery impedance, dc and ac voltage, ac current, frequency and battery temperature. The **BZA100** can obtain a Nyquist plot & Bode plot by measuring the impedance of the battery. This makes it easy to see changes in the R_s and R_{ct} values that correlate to the battery's state of health(SoH) and state of charge(SoC), allowing user to evaluate battery performance.


EIS data from **BZA100** can be analyzed with ZMAN impedance analysis software by automatic model searching and automatic fitting. Proper model library for user's batteries can be grouped to minimize the analysis time.

The user-friendly interface, compact design and rugged construction ensure optimal performance, test results and reliability. If you are familiar with java script coding, you can make your own software using web control.

BZA 100 Battery Impedance Analyser

Main Screen (Before the Experiment)

- Tab to see Bode, Nyquist or raw data list
- Real-time monitoring of current/voltage range, measured voltage value, and measured temperature value regardless of if a test is started. (data are not logged.)
- Displaying frequency, impedance, phase date of current measured point

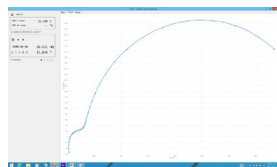

Parameter Input Box


Initial Frequency: 100000 Hz
 Final Frequency: 0.01 Hz
 Density: 10 Pts/decade
 Iteration: 1
 Current Range: 2 mA
 Max Initial Delay: 1.0 s
 Skip Cycle: 1
 Cycles: Auto

Parameter Input Box

- Galvanostatic EIS test
- Bias and amplitude current range is automatically set according to current setting range.
- Frequency range, number of data stored, number of repeated measurements, etc. can be set.



Readings During Measurement

- display a real-time Lissajous figure

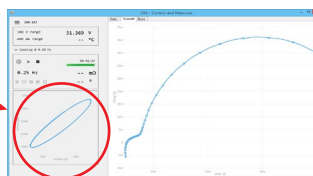


Main Screen After Measurement

- detailed information about data point (samples)

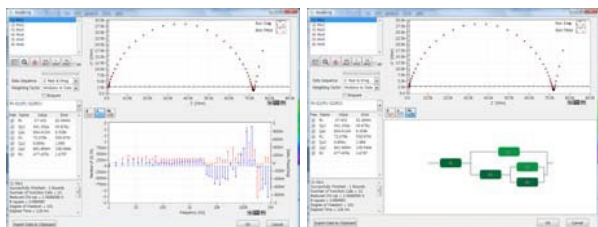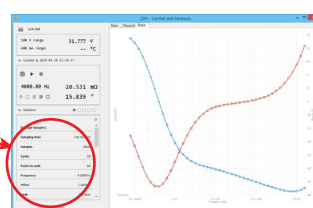

Fitting display

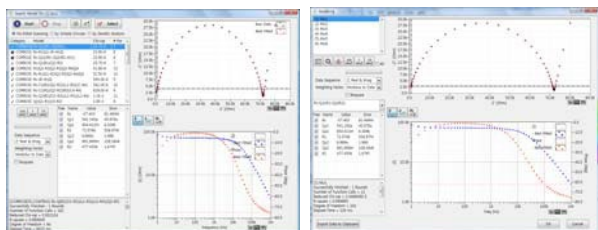

Automatic model searching



LEVM fitting

Standard cell cable with alligator clip



Option


- Large(Mid) size alligator clip
- Low impedance cable
- Cell cable modification
- Kelvin type small alligator clip cable (1 meter)
- Kelvin type medium alligator clip cable (1 meter)
- Pt100 Temperature sensor

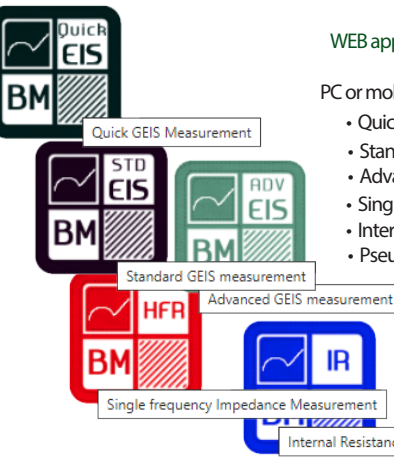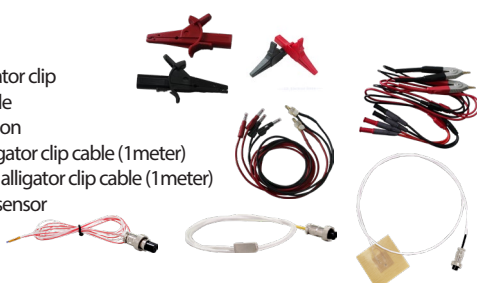

WEB applications

PC or mobile device control available

- Quick EIS measurement
- Standard EIS measurement
- Advanced EIS measurement
- Single frequency EIS measurement
- Internal Resistance measurement
- Pseudo Rs Rp measurement

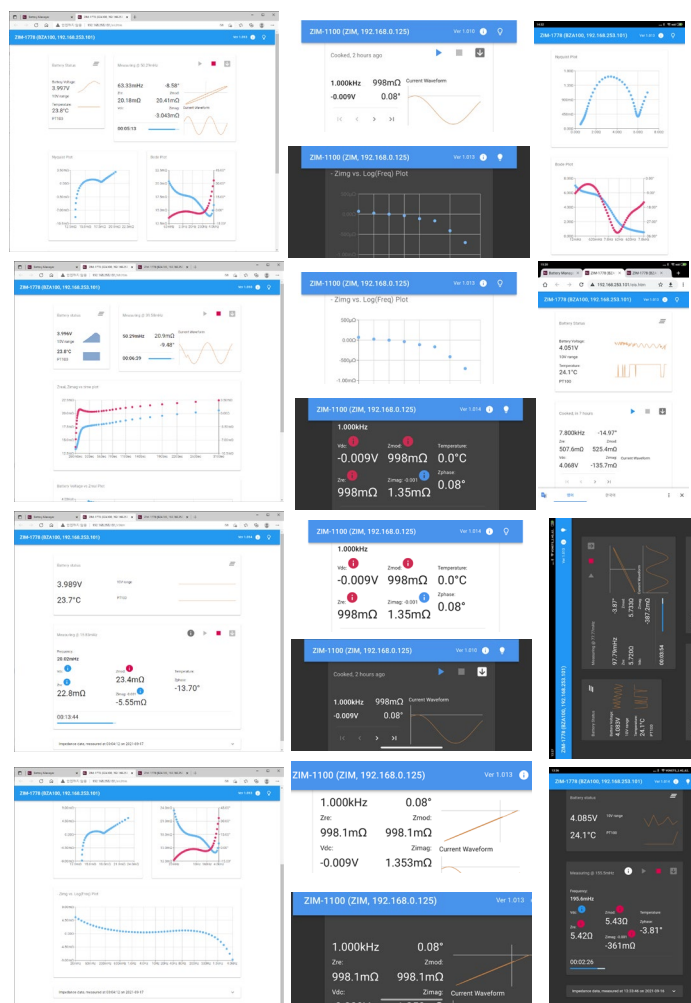


Note: For higher voltage than 100V battery, Select BZA1000 model for <1000V application.

High Current Cylindrical Battery Jig



Pouch Cell Jig

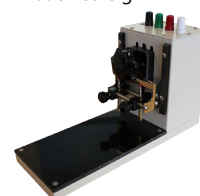

Cylindrical Battery or coin cell Jig



BZA 100 Battery Impedance Analyser

● Specifications

Impedance Measurement

Measurement range	500uΩ ~ 50Ω
Accuracy	±1% magnitude (1mΩ - 50Ω) ±1° phase
Frequency range	0.05Hz ~ 4kHz
Current amplitude (p-p)	400uA ~ 2A

DC Voltage Measurement

ADC resolution	24 bit
Input range	100V/10V (dual range)

AC Voltage Measurement

ADC resolution	24 bit
Input range	±250mV

AC Current Measurement

ADC resolution	24 bit
Current sensing Resistors	4ea (2A, 200mA, 20mA, 2mA)

Sinewave Generator

Frequency range	0.05Hz ~ 4KHz
Frequency accuracy	< 0.1%
Frequency resolution	0.01% or 5000 steps/decade
DAC resolution	10 bit
Output gain	2ea(X1, X0.2) total 8 current ranges (2A, 400mA, 200mA, 40mA, 20mA, 4mA, 2mA, 400uA)

Temperature Measurement

Input	RTD probe (PT100)
Accuracy	Max 1°C


Communication

Interface	LAN communication
-----------	-------------------

General

Size	160mm x 60mm x 180mm (WxHxD)
------	------------------------------

All specifications are subject to change without notice.

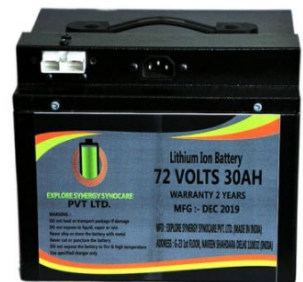

Designed by

ZIVE LAB
www.zivelab.com

WonATech

WonATech Co., Ltd.
7 Neunganmal 1-gil, Seocho-gu,
Seoul, 06801, Korea
Tel: +82-2-578-6516 Fax: +82-2-576-2635
e-mail: sales@wonatech.com
website: www.wonatech.com

Local Distributor